

Circulaire® Powder Containment



PCB 1800

**BRITISH
MANUFACTURED**



**FIVE YEAR
WARRANTY**



**GUARANTEED
PROTECTION**



**MAXIMUM
OPERATOR
SAFETY**




The Circulaire® Powder Containment Booth offers first class operator protection from both Powders and Particulate.

The unit is designed specifically for control of medium to heavy powders whilst dispensing from containment drums to a balance or smaller container. The high velocity rear extraction offers excellent containment by utilising H14 HEPA Filters, 99.997% efficient at 0.3µm and filtering air up to 1250m³/hr, achieving an air velocity of up to 0.7m/s.

Circulaire® Powder Weighing

The Circulaire® Powder Weighing Cabinet utilises Main & Outlet ULPA Filters to meet COSHH compliance. The unit is specifically designed for control of medium to heavy powders whilst dispensing to a balance or smaller container.

The Granite work surface ensures accurate weighing for sensitive applications and the removable and hinged baffles make for easy cleaning.

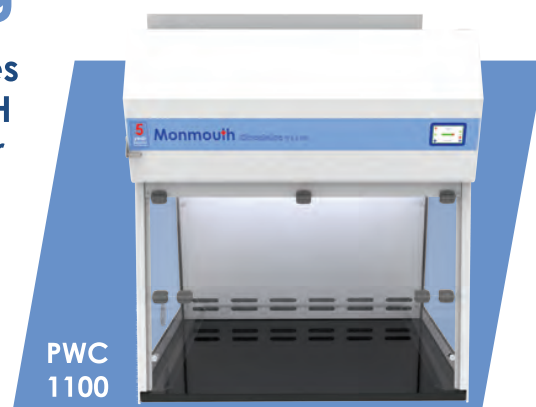

PWC
1100

Size Specification (W x D x H)

PCB 1800	PWC 1100
1800mm x 750mm x 2400mm	1100mm x 700mm x 1284mm